


# AUROCAUTERY

## Wet Field Bipolar Coagulator (Solid State)

Aurolab manufactures Aurocautery, which is a wet field bipolar coagulator. Mainly used in ophthalmic surgeries. It offers precise hemostasis of blood vessels like episcleral, intrascleral and intra ocular hemostasis and provides precise tissue coaptation like conjunctiva closure. Tissue may be coagulated under irrigation (wet - field) with saline or under bloody field.

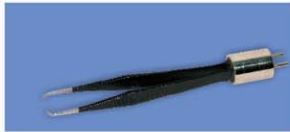

### SPECIFICATION

- Model : ACD-01
- Supply Voltage : 110v ac(on request)
- Radio frequency : 455 KHz
- Output RF Power : 7W approx
- Dimension : 140(W) x 50 (H) x 170(B) mm
- Weight : 1.5 Kg approx
- Casing : Powder coated M.S body

### FEATURES

- High frequency solid state circuitry
- Compact in size, elegant look and easy to use.
- Reusable forceps, erasers and cables
- Foot switch control,
- Audible indication
- Dial for variable power output
- Supplied with test certificate and Instructions for use leaflet.

Information Published in this catalogue are subject to change

 **aurolab**  
Excellence . . . in sight

Aurolab, No 1, Sivagangai Main Road, Veerapanjan, Madurai 625 020. India  
Phone: 91-452-2446100; Fax: 91-452-2446200  
E-mail: sales@aurolab.com; Web: www.aurolab.com

Issue No.3 - 07/08